



bagnoletto®



Allm**o**bility

BAGNOLETTO is the innovative **system that provides the pleasure of a real bath**. It ensures greater relief, both physically and mentally, safeguarding the primary goal of timely and accurate cleaning.

If you take care of a bedridden person, you surely know the importance of physical and psychological well-being and supporting the normal biological functions and the daily care of the body. Even when the patient is not self-sufficient, **it is essential to try to keep his hygienic habits unchanged**, supporting cleansing and hydration.



- ⊖ *The method currently used is sponging. This consists of washing the body with sponges (soaked in water or neutral detergent solutions), and drying the skin carefully. However, this method has many limits: sponges often cause irritation because of excessive rubbing (in these cases the skin is particularly sensitive). Furthermore, they do not always guarantee complete cleaning, especially in areas that are particularly difficult to reach (sores and interstices), and, if not properly sterilized they can be vehicles of infections (viruses and bacteria).*
- ⊖ *The other option, if the person's conditions allow it, could be the assisted bath and shower. A complete bath involves the use of a seat aid or of an electric lift, which must be installed on the bathtub, and the shower involves a wall seat or a wheelchair. Both methods are rather difficult to manage, they require more than one operator and can result in physical injuries or severe stress for the patient.*
- ⊕ **Allmobility Bagnoletto has been devised in order to overcome these problems.** It is the best alternative to traditional aids; it is **comfortable, practical, fast** and **ensures accurate cleaning**, reducing operator stress and discomfort for the disabled person.

 - ✓ **It can be used directly on the patient's bed, avoiding any incorrect handling of the patient**
 - ✓ **It offers all the benefits of a complete bath**
 - ✓ **It improves the quality of care and hygiene operations**
 - ✓ **It allows the use of running water in any conditions or context**
 - ✓ **No need of electricity**
 - ✓ **Suitable for any type of bed (even without sides)**
 - ✓ **Fast and easy to use**



BAGNOLETTO is a PVC sheet with stiffened side edges that are lifted and hooked together to form a container. It is positioned more or less in the same way as when changing the sheet with the patient on the bed. The operation is facilitated by the gripping straps.



Extend BAGNOLETTO

Position the person slowly on one side



Assembly

Place BAGNOLETTO partially under the person. Rotate the person and fully unfold BAGNOLETTO





Bathing

Use the pressure pump-tank and the hand shower. Give the pleasure of a complete bath



Dismantling

Empty the excess water. Remove the disposable sheet, if used, and disassemble BAGNOLETTO

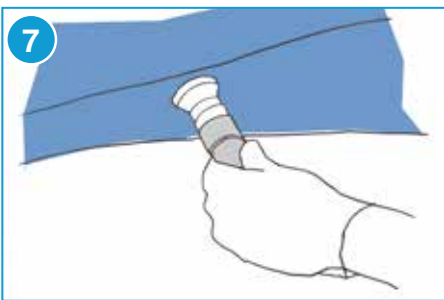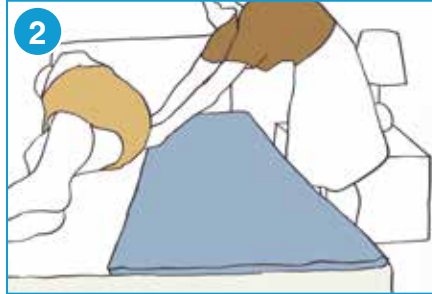

The system has several water drainage points and straight or curved corrugated hoses that facilitate this operation.

The water can be brought from a tap with an extension, or the manual pressure tank (accessory) can be used. For professional use in residential or rehabilitation facilities, it is suggested to employ disposable sheets that can be used for every single wash.



Watch the video on Youtube >



Shoulder carrying bag (included as standard). Useful for storing BAGNOLETTO

• **How to use BAGNOLETTO?**

Extend bagnoletto under your patient and use the pump-tank (which contains clean warm water) to shower the person. The pump tank is an accessory.

• **Can it be used with any type of bed?**

Yes, BAGNOLETTO is suitable for any type of bed (even without sides).

• **What are the dimensions of the BAGNOLETTO?**

85 x 190 x 20 cm

• **Which materials is BAGNOLETTO made of?**

BAGNOLETTO is made of polyester fabric and class 2C PVC.



• **How do I drain the waste water?**

After the bath, dirty water is drained into specific buckets (not provided), through the pipes provided.

• **How much water does the tank contain?**

The pump tank contains 10 litres.

• **Is it possible to add drugs to the water?**

Yes, it is possible to supply solutes and medical solvents.

• **Is it necessary to move the patient?**

No, BAGNOLETTO can be used directly on the person's bed and allows you to avoid any careless movement.

• **How many operators are necessary to use BAGNOLETTO?**

One operator is sufficient. The sequence is short and easy to follow, as shown in the video uploaded.

• **Is it possible to use BAGNOLETTO without electricity?**

Yes, BAGNOLETTO can be used in environments without electricity.

• **Is it necessary to sanitize the aid after use?**

Disposable sheets do not require any sanitizing, which makes BAGNOLETTO practical for healthcare facilities.



N.	ARTICLE CODE	DESCRIPTION
1	VQ BLBAGNOLE	BAGNOLETTO sheet 85x190x20 cm and draining pipes
2	VQ BLSERBAT	Manual pressure pump tank for BAGNOLETTO
3	VQ BLTELINI	Disposable sheets for BAGNOLETTO 10 pcs and 4 punches
4	VQ BLSOLOTELI	Disposable sheets for BAGNOLETTO 10 pcs
5	VQ BL RECIP	Folding bucket for water drainage Capacity: 11 lt



 bagnoletto®

