



EXACTA™

numberone **M**


EXACTA™ **BOLUS**
numberone **M**

EXACTA™ **UVA**
numberone **M**

EXACTA™ **DFO**
numberone **M**

EXACTA™
numberone **M**
MULTIRATE

EXACTA™
numberone **M**
MULTIRATE **BOLUS**

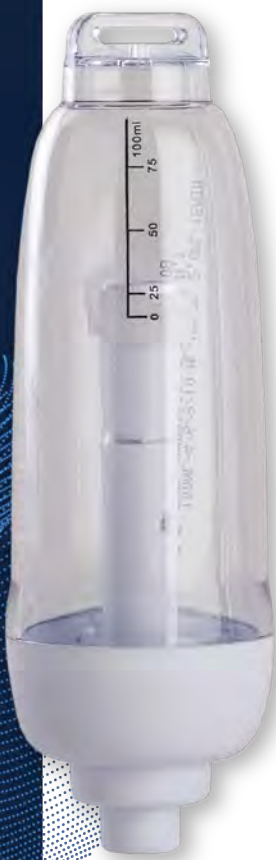

GAMASTECH
INNOVATIVE INFUSION SYSTEMS

The best choice for accurate, easy and safe continuous infusions



EXACTA™
number one **N1**

EXACTA™ DFO
number one **N1**
EXACTA™ UVA
number one **N1**



Designed to guarantee the absolute safety of intravenous, intra-arterial, epidural, subcutaneous, soft tissues and body cavities infusion therapies.

Recommended for the continuous administration of cancer, anaesthetic, iron chelation (DFO), antibiotic, and analgesic therapies and for the treatment of acute and chronic pain.

Suitable for home-care professionals.

In a few minutes it is easily filled up and ready for use.

- Hard, transparent reservoir with graduated scale
- Anti-kinking, triangular-shaped section infusion line in anti-UV material, DEHP FREE and LATEX FREE
- In-line filter with double filtering capacity of 1.2 microns (5 microns on request) for particles and antibacterial filtration and 0.02 micron air-venting
- Cap with hydrophobic filter to expel the air during filling phase and to prevent liquid leakage
- The UVA Model is shielded against UV rays and is equipped with a needleless valve to avoid leakage or drug aerosol during the filling phase.

Safe: All Exacta Number One pumps are tested one by one with a unique and exclusive Air Test method.

Available models:

Exacta N1, Exacta N1 Bolus, Exacta N1 UVA, Exacta N1 DFO

Available volumes: 60ml, 100ml, 150ml, 200ml, 275ml

Available flows: 0.5 ml/h to 500 ml/h

Exacta N1 is calibrated following the ISO 28620 standard conditions

Changed or altered conditions will affect flow accuracy

Viscosity:

Saline solution (NaCl)
or 5% glucose (D5W)

Temperature:

32C° (89.6 F)

Accuracy:

+/- 10% of the
nominal volume



EXACTA™ BOLUS
number one **N1**



**Bolus (PCA)
specifications:**

Available volumes:

60ml/100ml/275ml

Additional bolus:

0,5 ml/h, 1ml/h, 2ml/h

Min. release time:

8min/15min/
30min/60min

Kit:

- Travel bag
- Multiperforated Catheter
- Three-way stopcock

**Exacta N1 is also
available with the
following extras:**

- UV RAYS shielded
- Needleless valve
- 5 micron filter



Multiperforated Catheter



Needleless valve



Travel bag



EXACTA™
number one **MB**

MULTIRATE

MULTIRATE SPECIFICATIONS

VOLUMES: 60, 100, 150, 200, 275 ml

FLOW:

- 1-2-3-4 ml/h
- 2-3-4-5 ml/h
- 2-4-6-8 ml/h
- 4-6-8-10 ml/h
- 5-7-10-12 ml/h
- 6-8-10-12 ml/h

EXACTA™
number one **MB**

MULTIRATE

BOLUS

MULTIRATE WITH BOLUS SPECIFICATIONS

VOLUMES: 60, 100, 150, 200, 275 ml

FLOW:

- 1-2-3-4 ml/h
- 2-3-4-5 ml/h
- 2-4-6-8 ml/h
- 4-6-8-10 ml/h
- 5-7-10-12 ml/h
- 6-8-10-12 ml/h

LOCKOUT TIME:

8/15/30/60 min.

SAFETY

Flow Rates can be changed only by clinicians thanks to its special key.

DISCRETE

Designed to take up as little room as possible.

EASY TO USE

It's easy to change the flow rate during the infusion, thanks to its special key.

VERSABLE

More than 600 flow rates to meet the needs of the most important and common therapies.



Via Giovanni Verga, 27
95030 Sant'Agata Li Battiati (CT) - Italy
Tel. +39 0952938579 - Fax +39 0955183766
www.gamastech.com - info@gamastech.com

CE 0051

Distributed:

Rev. 11/2019