



What does dual-principle loop mean?

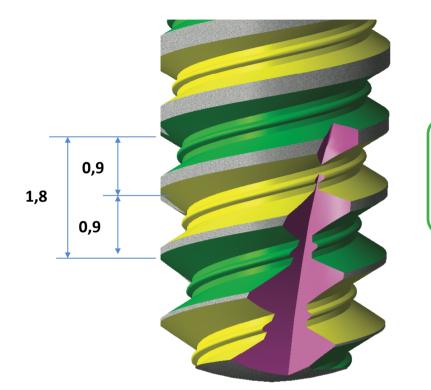

Easy Slim has a pitch of 0.9 mm
The system makes a turn and goes in 1.8 mm

The **Easy Slim** system has a double spiral

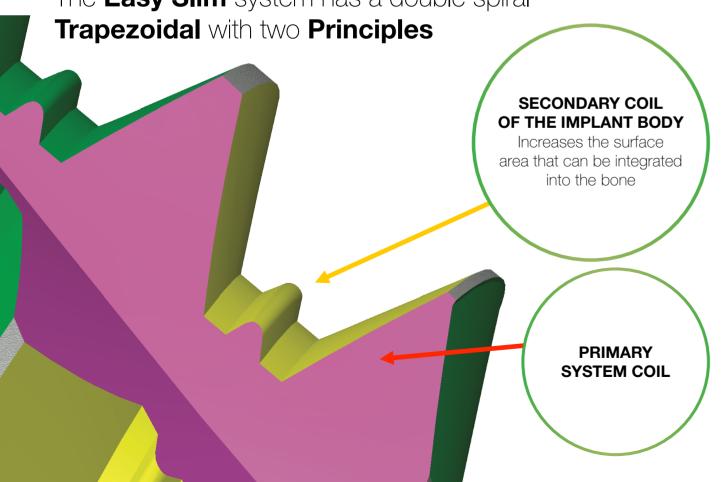



Ornaghi s.r.l.

Viale Marche, 70 - 20159 Milano (MI) Phone +39 0392142070 www.ornaghimedical.com












EasySlim


Implantline





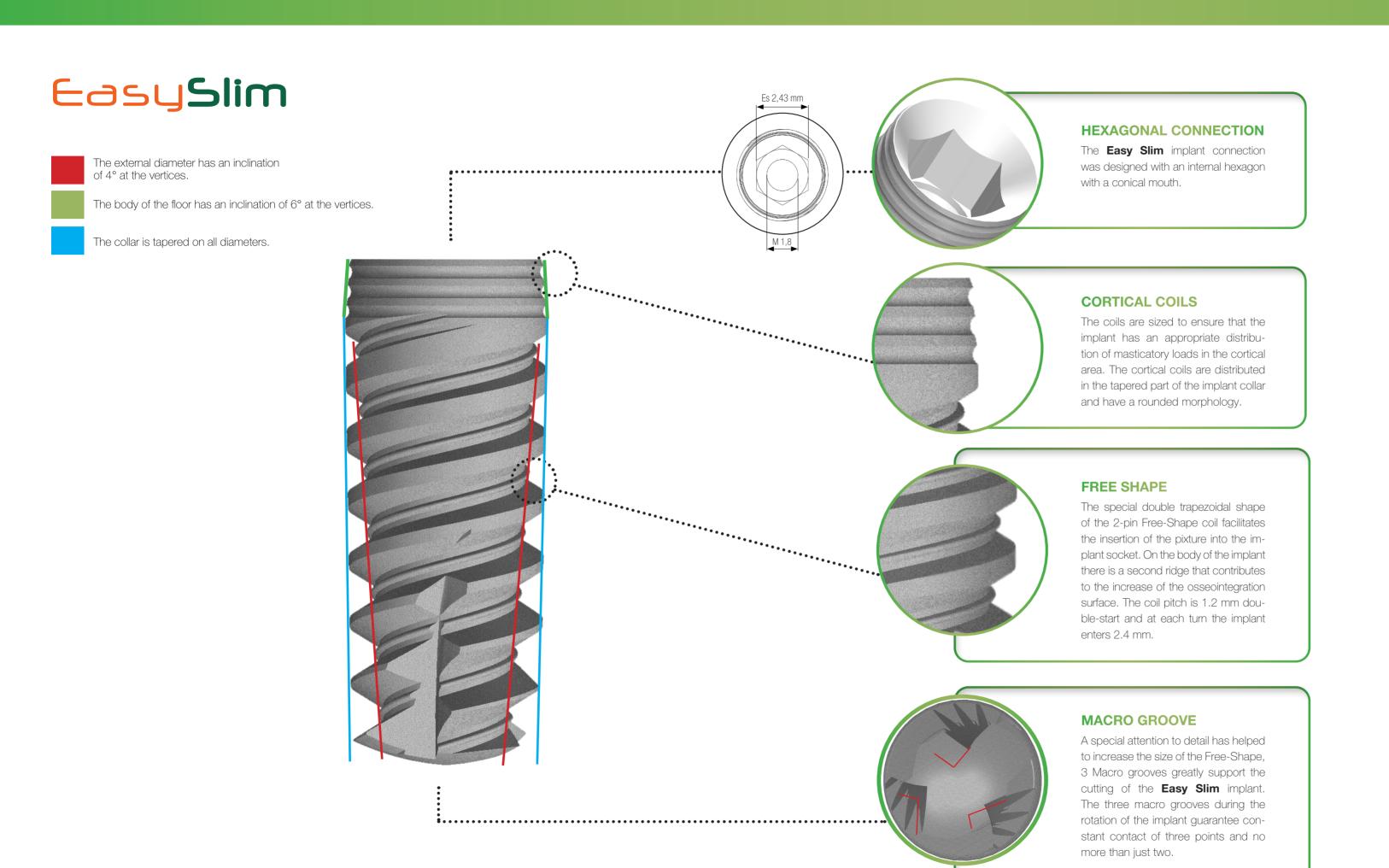


The system has a double principle coil

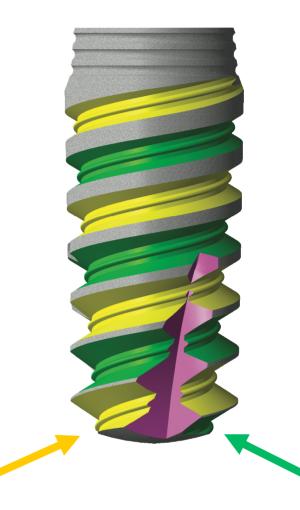

Double Principle means that the system has two threads and each starts at 180° from the other

In the image we see
the first is yellow and the second
green coil which is 180° in diameter
(opposite side)

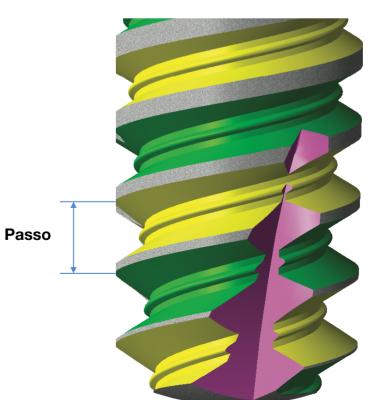

Means that the system at each turn screws twice the value of the *Passo della Spora*.

The wheelbase of the implant is the distance between two points taken in the same position on two ridges. **Easy Slim** has a 0.9 mm pitch